

1.COLD TREATEDWATER INTAKE OF SUPPLIED AT TEMPERATURE OF 25°C, CONSUMPTION PER 8~10L/CYC, HARDNESS OF WATER≦0.2°fH(0.02mmol/L),pH 5~7AND WITH CONDUCTIVITY ≤ 15 µ S/cm

E ACH TIME THE STEAM GENERATOR IS EMPTIED FOR CLEANING OR MAINTENANCE, A FURTHER 15 LITER OF WATER ARE NEEDED TO FILL IT UP AGAIN. 2. DRAINAGE PIPE Ø≥1" METALLIC PIPE, CORROSION RESISTANCE 2"G. AND A

SELECTRIC INTAKE ,12kW, WIRE 4 m LONG. (GENERAL POWER SUPPLY TO THE STERILIZER) TEMPERATURE OF 100°C ABOVE.

<u></u>		<u> </u>
440VA	1	15.7A
415VAC		16.7A
380VAC		18.2A
230VAC 380VAC 415VAC 440VAC	52.2A	30A
	3₩	41/
	1ϕ	34

A. AIR INTAKE FREE OF IMPURITIES #1/4" FILTERED AND REGULATED AT A PRESSURE OF 6-7 bar, FLOW OF 0.0083 Nm³/min, APROX. CYCLE CONSUMPTION

0.005 Nm³/min, PLACED THE WALL WITH BUILT-IN AND ACCESSIBLE KEY.

THIS INTAKE WILL BE PROVIDED ONLY IF THE CUSTOMER HAS AIR NETWORK

TO MODIFY THIS DRAWING ACCORDING TO

NOTE

FHIS COMPANY RESERVES THE RIGHT TO MODIFY TI THE TECHNICAL DEVELOPMENT OF THE EQUIPMENT.

IN ACCORDANCE WITH ESTABLISHED APPLICABLE NORMS AND REGULATIONS. THE COMPAY DECLINES ALL RESPONSABILITY FOR FAILURE TO COMPLY WITH SUCH LAWS AND REGULATIONS.

/I/WARNING

NOTE: ONCE THERE EXISTS REVERSING OF FAN BLADES THAT WILL CAUSED NO VACUUM VALUE, . TURN OFF THE MAIN POWER AND THEN SWITCH ANY TWO LINES OF THE THREE-PHASES NOTE:POWER INPUT MUST BE INSTALLED ACCORDING TO LOCAL APPLICABLE NORMS. $\sqrt{}$

EXCEPT GROUND LINE. YOU MAY HAVE TO REPEAT THIS SWITCHING STEPS FOR TWO TIMES TO MATCH

WITH THE CORRECT ROTATING DIRECTION AS SHOWN ON THE PUMP. 4. CEILING MOUNTED AIR EXTRACTION WITH GRILLE DRAW 300 Kcal/h (LOADING)

DRAINAGE PIPE AND EXHAUST CONNECTING PIPE MUST KEEP THE STRAIGHT, UNLOADING AREA)

SHALL NOTBE BEN'

B. THE EQUIPMENT SHALL BE ADJUST TO A ANGLE OF 1°-2° WITH REAR SITE

HIGHER THAN FRONT SITE FOR THE PROPER OPERATION OF DRY AND DRAINAGE

MPORTANT: ENVIRONMENT CONDITIONS IN THE STERILE AREA:

TO KEP THESE CONDITIONS IN THE STERILE AREA AND DEPENDING ON THE KCAL LOST BY THE EQUIPMENT INDICATED IN THE ITEMS 4, WE SUGGEST THAT THE EXTRACTOR AND THE DIAMETER OF THE EXTRACTION TUBE ARE CALCULATED THE STERILE AREA TEMPERATURE MUST BE ≤ 40°C ;HUMIDITY ≤ 85%RH AND MADE BY A SPECIALIZED COMPANY

SELECTRIC INTAKE BOX

FALSE CELILING LINE

800

BAIR INTAKE


DEXHAUST FAN MOTOR

CELILING LINE


- THE ELECTRIC WIRES TO THE STERILIZER SHOULD BE DRIVEN BY AN HERMETIC CONDUIT
- THE INDICATED AIR EXTRACTION IS ONLY TO DISSIPATE THE HEAT OF THE MACHINE AND MUST BE INDEPENDENT OF THE CLIMATE CONTROL SYSTEM OF THE
- CHAMBER DIMENSIONS: Ø610x1200 mm
- DEVICE OUTSIDE DIMENSIONS: W800mmxH1720mmxD1500mm ALL THE OUTLETS MUST BE INSTALLED ACCORDING 70 LOCAL APPLICABLE NORMS IF THERE IS NO FALSE CEILING, THE AIR EXTRACTION, MUST BE INSTALLED IN THE WALL, AND HAVE THE SAME CAPACITY AS THOSE OF THE FALSE CEILING, SO THE AIR CAN BE CARRIED OUTSIDE OF THE BUILDING. IT ALSO NEEDS TO BE EQUIPPED
- WITH A GRILLE. THE INSTALLATION UP TO THE MACHINE AS WELL AS THE SUPPLY OF ACCESSIBLE KEY ACCORDING TO THE SPECIFICATIONS INDICATED IN THE DRAWING, ARE
- IF THERE IS ANY KIND OF PROBLEM DURING INSTALLATION FOR THE VARIOUS INTAKES NOT INDICATED IN THE DRAWING, PLEASE CONTACT OUR TECHNICAL SERVICE DEPT. FOR FURTHER ASSISTANCE.

CHARGEABLE TO THE PURCHASERS.

- THE MEASUREMENTS OF THIS DRAWING ARE CONSIDERING FINISHED WALLS.
- MIN. SPACE FOR SERVICE. THE SPACES OF MAINTENANCE MUST BE SUFFICIENTLY ILLUMINATED.

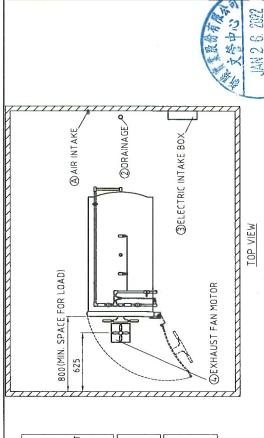

SELECTRIC INTAKE BOX MIN. SPACE FOR SERVICE **BAIR INTAKE** 800 Ξ̈́ CELILING LINE SE CELILING 2979(STERILIZE SPACE) *EXHAUST FAN MOTOR* 1500 800 625

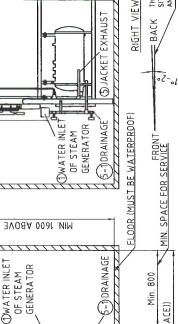

.... 00

1720

0081

FINISHED WALL > 2000mm

(PIPE-END HIGH ≤ 100mm

ODRAINAGE PIPE (METAL Ø≥ 1") THE EQUIPMENT SHALL BE ADJUST TO A ANGLE OF 1*-2*WITH REAR SITE HIGHER THAN FRONT SITE FOR THE PROPER OPERATION OF DRY AND DRAINAGE.

BACK

FRONT VIEW

SECTION(1:5)

®LEVEL ÁDJUSTER

parts should comply with RoHS.

Note

JAN.26.2022

文節中心 案股份有

2400((STERILIZE SPACE))

B SECTION

(5-2) CHAMBER EXHAUST 800

Ξ̈́

PLAN MARKA SA-600A INSTALLATION
STATES SA-600A-HE-001-03 110 2505 STURDY INDUSTRIAL なななる。 20211221 Denny DR WN 参 20120711 Joe LD, 未標示之圖角為KI導角為CI + 限公 ANGLE MODEL NO. 许表 機 编 SA-600A 業股份有 501~X ±2 201~500±1 1:40 SCALE 比例 新駿實 OLERANCE未標註者一般公差表 101~200±0.8 51~100±0.6 JRDY 1~50 ±0.4 SURFACE FINISHING 表面處理 0~10 ±0.2

1020

73-058-06