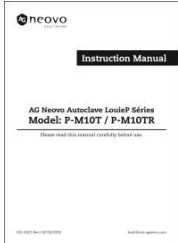
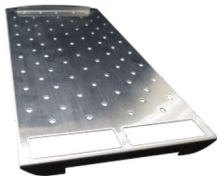







Packing List

Packed Items	P-M10TR
 <p>Instruction Manual (1)</p>	✓
 <p>Heater Cover (1)</p>	✓
 <p>Tray Rack (1)</p>	✓
 <p>Trays (3)</p>	✓
 <p>Tray Holder (1)</p>	✓
 <p>USB Flash Drive (1)</p>	✓
 <p>Cleanser for Chamber (4)</p>	✓

Ambient Temperature Range	41 °F to 104 °F (5 °C to 40 °C).
Relative Humidity	80% RH @ 88 °F (31 °C) to 50% RH @ 104 °F (40 °C). (Pollution Degree 2, in accordance to IEC 60664)
Operating Attitude	< 3,280.84 ft. (1,000 m)
Device approved for Indoor Use Only	
Dimension (D x W x H)	25.8 in. x 21 in. x 17.4 in. (65.5 cm x 53.3 cm x 44.2 cm)
Weight	Net weight : 127.9 lbs (58 kg) Gross weight : 134.5 lbs (61 kg)
Voltage Fluctuations	±10%
Electrical Requirement	12A / 120VAC 60Hz / 1400W

Note: Product specifications and accessories may vary and may be subject to change without prior notice.

Location Requirements

1. It will require at least two (2) or more people to carry the sterilizer to avoid dropping it off by mistake.
2. While installation, please make sure that the bearing capacity of installation table is enough to sustain the weight of the sterilizer.
3. Position the sterilizer on a stable bench, laboratory bench or work surface, ensuring at least **4 inches (10 cm)** clearance between the wall or other pieces of equipment and the sides of the unit for free circulation of air.
4. Make sure the sterilizer door can be opened freely after installation.
5. Do not install or operate the sterilizer in areas where flammable items or volatile substances are used or stored. An explosion could occur, causing personal injury. An installation site with good air circulation is required.
6. A separate (dedicated) socket is required for the sterilizer. Make sure the socket is earthed and can offer the power source with **12A / 120VAC 60Hz / 1400W**.

