

### WARNING

THIS DRAWING GIVES THE BASIC INSTRUCTIONS FOR PROVIDING THE SUPPLY SERVICES AND ENVIRONMENT THAT GUARANTEES THE PROPER FUNCTIONING OF THE PRODUCTS FROM THIS COMPANY. THE PURCHASER, AND THE PURCHASER'S INSTALLERS SHALL PROVIDE THE SUPPLY LINES NECESSARY FOR EACH PRODUCT IN ACCORDANCE WITH ESTABLISHED APPLICABLE NORMS AND REGULATIONS. THE COMPANY DECLINES ALL RESPONSIBILITY FOR FAILURE TO COMPLY WITH SUCH LAWS AND REGULATIONS.

### NOTE

THIS COMPANY RESERVES THE RIGHT TO MODIFY THIS DRAWING ACCORDING TO THE TECHNICAL DEVELOPMENT OF THE EQUIPMENT.

This equipment has been designed for use in accordance with the International EMC (Electromagnetic Compatibility) Standards. In view of different environments, please follow the instructions given below to eliminate interference, if necessary.

- Move the equipment or rotate its direction;
- Enlarge the space between the equipment and other machines;
- Put the plug into other outlets;
- Please consult with the local distributor or qualified electrician.

### Installation :

1. It will require at least two (2) or more people to carry the sterilizer to avoid dropping it off by mistake.
2. While installation; please make sure that the bearing capacity of installation table is enough to carry the sterilizer. For the weight information of the sterilizer, please refer to "Specifications".
3. Located position the sterilizer on a stable bench or work surface.
4. Position the sterilizer on a stable laboratory bench or work surface, ensuring at least 10 cm clearance between the wall or other pieces of equipment and the sides of the unit for free circulation of air.
5. Make sure that the door can be opened freely after installation.
6. Do not install or operate the sterilizer in areas where flammable items or volatile substances are used or stored. An explosion could occur, causing personal injury. An installation site with good air circulation is required.
7. ELECTRIC INTAKE : 1  $\phi$  / 220V/3000W/13A WIRE 4 m LONG.  
(GENERAL POWER SUPPLY TO THE STERILIZER).

### ENVIRONMENT CONDITIONS IN THE STERILE AREA:


- Indoor use;
- Under 1,000m (altitude);
- Temperature 5°C to 40°C;
- Relative Humidity 80%RH@31°C to Relative Humidity 50%RH@40°C;
- Voltage fluctuation  $\pm 10\%$ ;
- Transient overvoltages category III;
- (Pollution degree 2)

### Important Safety Instructions :


1. Connect the external water supplier to the water inlet with pressure greater than 2 bar.
2. Refer to the table for the requirements of water quality.
3. Connect the drainage pipe to exhaust the used water according to the local authority's law.
4. Do not obstruct the drainage system.
5. Consult technicians for the correct connections.

### ATTENTION:


- CHAMBER DIMENSIONS:  $\phi 300$  mm x 450 mm.
- DEVICE OUTSIDE DIMENSIONS: 584(W) mm x 480(H) mm x 670(D) mm.
- ALL THE OUTLETS MUST BE INSTALLED ACCORDING TO LOCAL APPLICABLE NORMS.
- THE INSTALLATION UP TO THE MACHINE AS WELL AS THE SUPPLY OF ACCESSIBLE KEY ACCORDING TO THE SPECIFICATIONS INDICATED IN THE DRAWING, ARE CHARGEABLE TO THE PURCHASERS.
- THE MEASUREMENTS OF THIS DRAWING ARE CONSIDERING FINISHED WALLS.
- THE SPACES OF MAINTENANCE MUST BE SUFFICIENTLY ILLUMINATED.

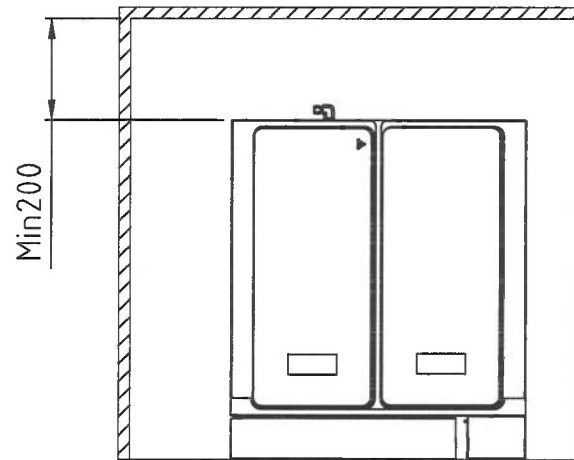

TOP VIEW

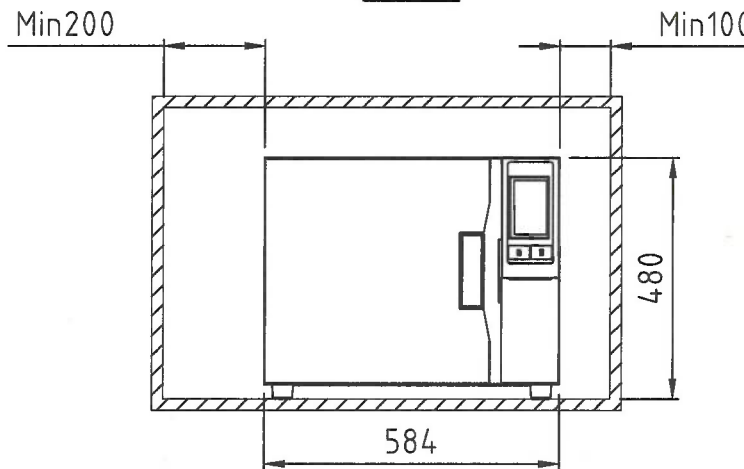

FRONT VIEW

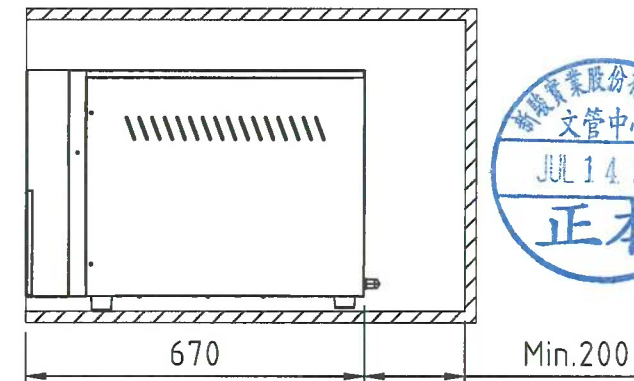

RIGHT VIEW

### NOTE

Installer Qualifications: The install personnel having proper qualification and holding technical documentation (including Technician Manual) and adequate information are authorized to install the apparatus.

Table : Suggested maximum limits of contaminants in and specification for water for steam sterilization

	feed water
Evaporate residue	$\leq 10$ mg/l
Silicium oxide ,SiO <sub>2</sub>	$\leq 1$ mg/l
Iron	$\leq 0.2$ mg/l
Cadmium	$\leq 0.005$ mg/l
Lead	$\leq 0.05$ mg/l
Rest of heavy metals ,excluding iron, cadmium,lead	$\leq 0.1$ mg/l
Chloride	$\leq 2$ mg/l
Phosphate	$\leq 0.5$ mg/l
Conductivity(at20°C)	$\leq 15$ $\mu$ s/cm
pH value	5 to 7.5
Appearance	colourless, clean, without sediment
Hardness	$\leq 0.02$ mmol/l

NOTE The use of water for steam generation with contaminants at levels exceeding those given in this Table can greatly shorten the working life of a sterilizer and can invalidate the manufacture's warranty of guarantee.



## Sturdy 新駿實業股份有限公司 STURDY INDUSTRIAL CO., LTD.

TOLERANCE	未標註者一般公差表	未標示之圓角為R1導角為C1	DESIGNED	DRAWN	VERIFIED	APPROVED
0~10 $\pm 0.2$	51~100 $\pm 0.6$	201~500 $\pm 1$	設計: Mia	繪圖: Mia	校核: 侯雲聖	核准: Rohf
11~50 $\pm 0.4$	101~200 $\pm 0.8$	501~X $\pm 2$	20220714	20220714	20220714	2022.7.14
MATERIAL	SCALE	VIEWPOINT	NAME	NAME	NAME	NAME
材質	比例 1:15	視角	圖名 PRIME30B INSTALLATION PLAN	圖名 PRIME30B INSTALLATION PLAN	圖名 PRIME30B INSTALLATION PLAN	圖名 PRIME30B INSTALLATION PLAN
SURFACE FINISHING	MODEL NO.	UNIT	DWG. NO.	DWG. NO.	DWG. NO.	DWG. NO.
表面處理	代表機種 PRIME30B	單位 mm	圖號 PRIME30B-HE-001-01	圖號 PRIME30B-HE-001-01	圖號 PRIME30B-HE-001-01	圖號 PRIME30B-HE-001-01